


MEDICINE AND AESTHETICS

Motus AX



MOTUS AX

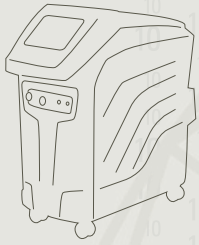

The New Era of Hair Removal

The missing solution!

- Hair Removal
- Benign Pigmented Lesions



The Code of Excellence



THE MISSING SOLUTION IN THE HAIR REMOVAL FIELD

The hair removal market is extremely complex. The numerous brands presenting their products generate fierce competition. Every clinic tries to offer the best care for their patients, but it is troublesome to discern the difference among these systems.

All this is due to the fact that no definitive solution has yet been found for all patients and on all skin and hair types.

Until now!

DEKA is proud to introduce the new **Motus AX** Alexandrite laser system that has been specially designed and developed with Moveo technology that will revolutionise the laser hair removal market.

"I have been using the Alexandrite laser for years in my clinics because, in my experience, I consider it the most effective device for hair removal, especially when I have to remove fine and blond hair. No other system or wavelength can achieve satisfying results for my patients and myself in the same way. I do not deny, however, that in the past, other systems were simpler or more economical to use, but at the expense of a highly satisfactory final result. Since DEKA introduced the Moveo technology, and now even more so with the Motus AX, I have found the perfect product for my clinics, which satisfies my every need: a surprisingly simple method to learn and to perform; excellent results like never before with no more problems treating darker skin types or in the presence of light and fine hair; no complaints or discomfort by patients thanks to a completely painless treatment which has no side effects."

Prof. Nicola Zerbinati, M.D.

Department of Dermatology
University of Insubria - Varese, Italy



Motus AX



A NEW ALEXANDRITE LASER CONCEPT

Motus AX is the new system that DEKA has developed to complete the Moveo technology. A system without rivalries, with all the benefits of Alexandrite laser without the disadvantages:

- The only wavelength which is effective at treating fine and light hair. In fact, no other type of laser is able to combine an **efficient melanin absorption** (at 755 nm), a **high peak power** (more than 6,000 W) and the possibility of selecting a short pulse length (up to 2 ms).
- The Moveo handpiece allows unwanted hair to be removed by means of a new method that is easier, faster, safer and more effective (painless treatment, also suitable for darker skin types).
- An innovative technology that structurally simplifies the heart of the laser system, reducing costs and maintenance problems.
- A compact, portable system.
- A surprisingly low cost!

DEKA has changed the whole concept of Alexandrite laser with the **Motus AX**.



PLUS

755 nm	Wavelength universally recognised as the most effective for hair removal.
Moveo Technology	The new frontier for hair removal.
20 mm	One of the largest spot sizes for Alexandrite.
10	Choose up to ten handpieces on board for the laser sources (customized selection) with an automatic recognition system.
Cooling	Integrated contact skin cooling system.
Over 200	Integrated protocols developed for Moveo Technology and standard handpieces.
16 A	The Motus AX's electrical supply is limited to 16 A. The system can be plugged into any common electrical outlet in a clinic or doctor's office.
Simple & User Friendly	The Motus AX has been entirely designed by DEKA R&D team, which integrates elements of simplicity and maneuverability: compact size, intuitive software, lightweight and ergonomic handpiece with contact cooling, and the entire Moveo method itself.
Great Benefits	A rapid return on investment and a total patient satisfaction will increase your business!





MOVEO TECHNOLOGY: SIMPLY EXTRAORDINARY!




DEKA came out with a revolutionary technology in the field of laser hair removal. Moveo is not just a new product. A new idea based on 4 fundamental concepts: **effectiveness, speed, safety and simplicity.**

The Effectiveness of Technological Innovation

The energy emitted by a laser system with wavelength visible and near infrared radiation is largely reflected by the skin. The graph below shows the reflectivity of the skin at various wavelengths and melanin content. For instance, the 755 nm Alexandrite laser wavelength demonstrates the percentage of reflection with melanin content of 8% (corresponds to phototype II) is over 50%. This is a significant portion of energy lost during the treatment, which cannot be used for therapeutic purposes. This is due to the high variation of reflectance between the skin and the surrounding air.



The Moveo handpiece has a built-in skin contact cooling system. The sapphire tip is cooled to 15°C for maximum patient comfort.

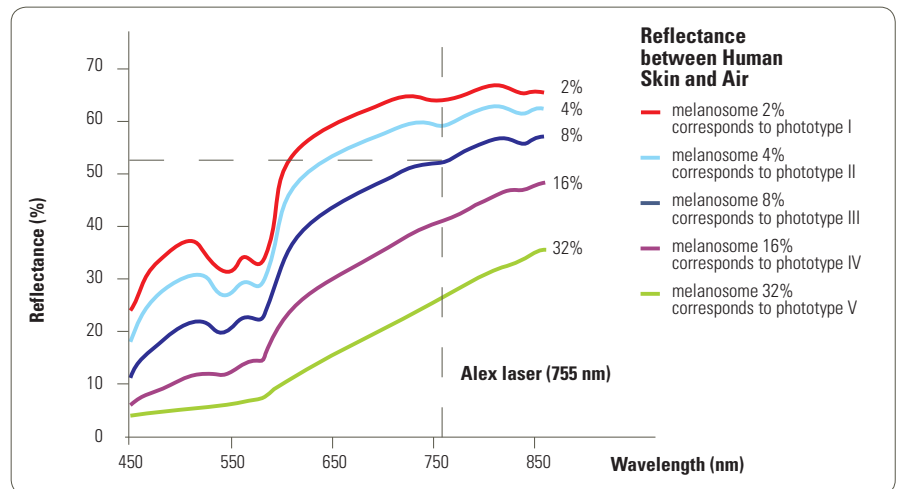

The Moveo handpiece was designed to solve this problem. DEKA has developed a technology that optimizes the laser-skin coupling by doubling the transmission of energy, with **an energy efficiency that has never been seen until now!** The sapphire tip that comes in contact with the skin decreases the variation in the reflection index by *reducing the reflected energy loss.*

All the Advantages but without Limitations!

DEKA has overcome limitations and difficulties of the traditional Alexandrite for hair removal while retaining all the advantages. The Moveo Technology uses simple, even movements over specific areas, **quickly** removing unwanted hair due to its **20 mm sapphire tip size.** The repeated passes over small areas, cause gradual heating of the vital parts of the hair leading to its destruction in a way that is totally **painless for the patient and side-effects free such as burning or discolouring.**

▶ ONE METHOD FOR ALL: EASY TO USE. EASY TO DELEGATE.




From an operational viewpoint the procedure is very simple:

- The handpiece is automatically recognized by the system that sets the Moveo mode.
- The **simplified software** and integrated database permit rapid choice of treatment parameters for any phototype, area and hair type.
- An accumulator displays the energy delivered during the treatment and an acoustic signal warns the user when the desired end-point has been reached. At this point you can treat another skin area.
- With the magnitudes involved (fluence, frequency and energy) the treatment time for a 10 cm x 10 cm area is approximatively 25 seconds.

The Moveo technology is **quick, easy to learn** and can be put **immediately into action** by different categories of staff (doctors, medical assistants, nurses).

▶ CLINICAL CASES *[Courtesy of P. Bonan, M.D. and M. Troiano, M.D. – Florence, Italy]*

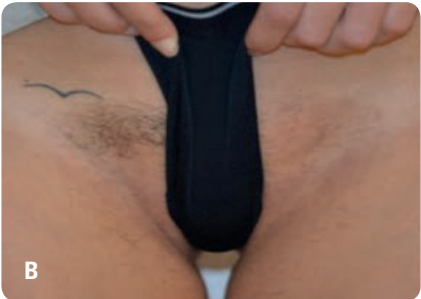

Bikini line treated with **Moveo Technology**. The patient had a small black tattoo on her right side. For precautionary reasons the doctor left untreated the area close to the tattoo (2-3 cm) during the first session. After a month and a half from the only treatment administrated, it is easy to observe the treated area (hairless), comparing to the hair regrowth adjacent the tattoo area.



Male back hair removal treatment. Only left side was treated once with Motus AX laser system using **Moveo** handpiece. The picture shows on the left side (A) the condition 5 months after one session, while the right side (B) shows the initial patient condition without any treatment.

TECHNICAL DATA

Motus AX

Laser Type	Alexandrite
Wavelength	755 nm
Energy per Pulse	31 J (Max.)
Laser Power	62 W (Max.)
Number of Pulses	1 or 2
Pulse Duration	From 2 to 50 ms (single pulse) From 4 to 80 ms (twin pulse)
Total Pulse Length	From 2 to 300 ms
Repetition Rate	From 1 to 10 Hz
Spot Sizes (optional)	5 mm, 7 mm, 10 mm, 12 mm, 14 mm, 15 mm, 16 mm, 18 mm, 20 mm; moveo handpiece
Fluence	Up to 160 J/cm ²
Aiming Beam	Green Diode Laser @ 532 nm < 4 mW
Emission Control	Fingerswitch or footswitch
Laser Handpiece Cooling	Integrated Skin Cooler
Dimensions	95 cm (H) x 51 cm (W) x 83 cm (D)
Weight	90 kg
Electrical Requirements	230 Vac - single phase - 50/60 Hz - 16 A

CAUTION

Visible and invisible laser radiation.
Avoid eye or skin exposure to direct or scattered radiation.
Class 4 laser product.



This brochure is not intended for the market of USA.

Hair Removal - Benign Pigmented Lesions

CE
0459



The Code of Excellence

Follow us on



www.dekalaser.com

Motus AX

Dealer stamp



DEKA M.E.L.A. s.r.l.
Via Baldanzese, 17 - 50041 Calenzano (FI) - Italy
Tel. +39 055 8874942 - Fax +39 055 8832884

DEKA The Code of Excellence

A spin-off of the EI.En. Group, DEKA is a world-class leader in the design and manufacture of lasers and light sources for applications in the medical field. DEKA markets its devices in more than 80 countries throughout an extensive network of international distributors as well as direct offices in Italy, France, Germany, Japan and USA. Excellence is the hallmark of DEKA's experience and recognition garnered in the sphere of R&D in over thirty years of activity. Quality, innovation and technological excellence place DEKA and its products in a unique and distinguished position in the global arena. DEKA manufactures laser devices in compliance with the specifications of Directive 93/42/EEC and its quality assurance system, certified by G.I.S. in accordance with the ISO 9001 and ISO 13485 standards.