

# EBS RHAPSODY

## 3000W

### 755nm+808nm+1064nm

Super High Power Diode Laser Hair Removal Machine




Gold Standard working energy- 10J + 5ms

Super short pulse duration

Super long life span

Super high power

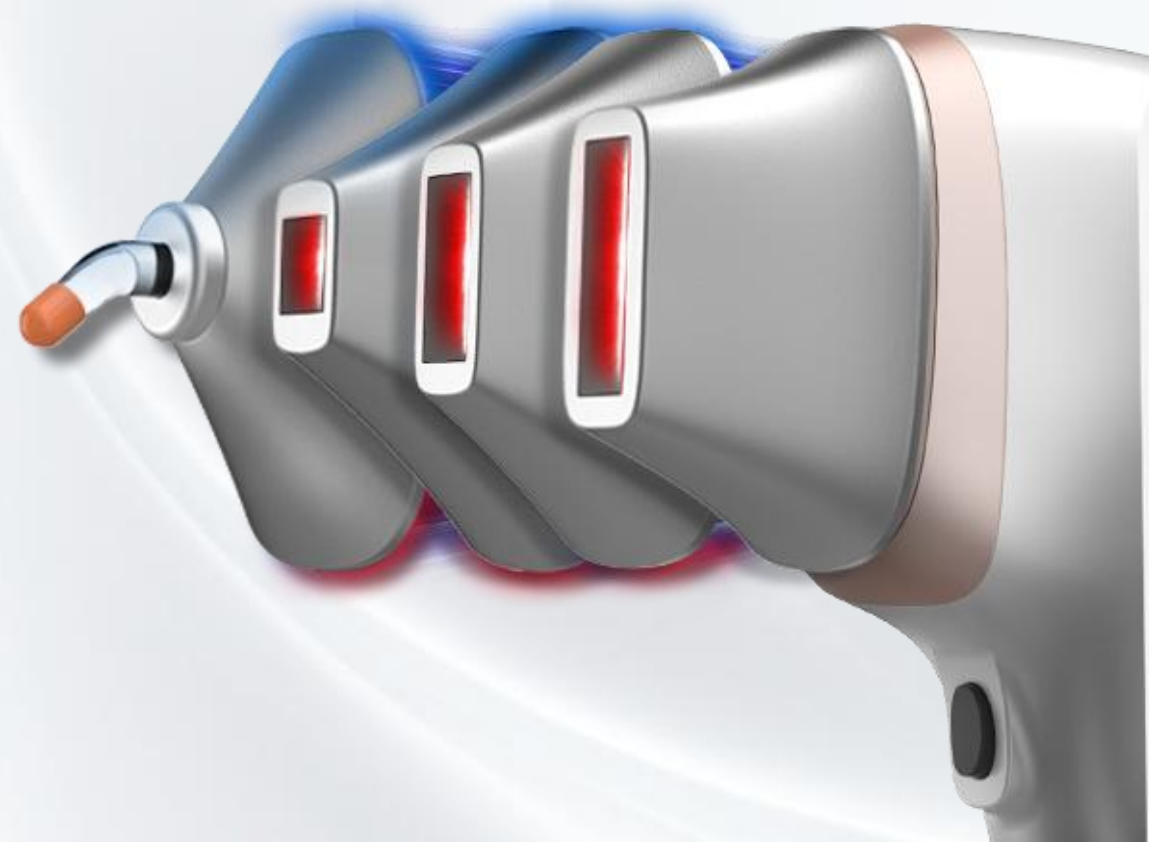

# Product Display



Handpiece With Touch Screen

12inch Big Screen

Holder

Emergency Button

Switch Button


High-End Hollow Silent Wheel




# Treatment Head Advantage

Different spot sizes head working on one handle. More Choices for different body area.

**Standard**

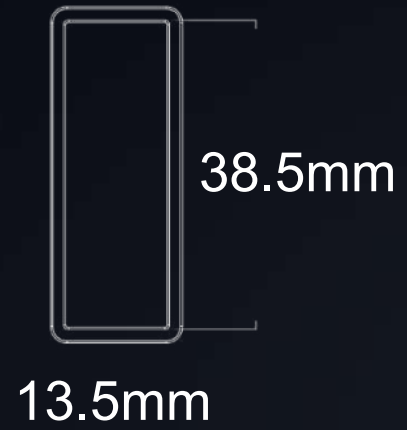

**Standard**

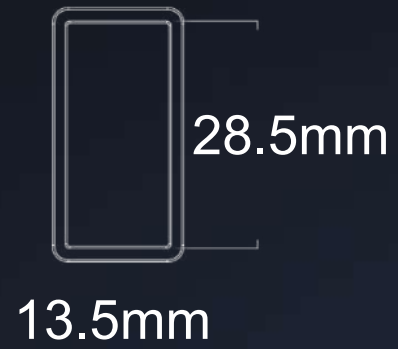

**Standard**



**Standard**



# Treatment Site

**Face:**  
12.5\*14.5mm or Round 8.25mm

**Forehead:**  
12.5\*14.5mm

**Nose Hair:** Round 8.25mm  
**Mustache:** 12.5\*14.5mm

**Arm:**  
13.5\*28.5mm


**Armpit:**  
12.5\*14.5mm

**Chest:**  
13.5\*38.5mm

**Back:**  
13.5\*38.5mm

**Bikini Line:**  
12.5\*14.5mm

**Leg:**  
13.5\*28.5mm

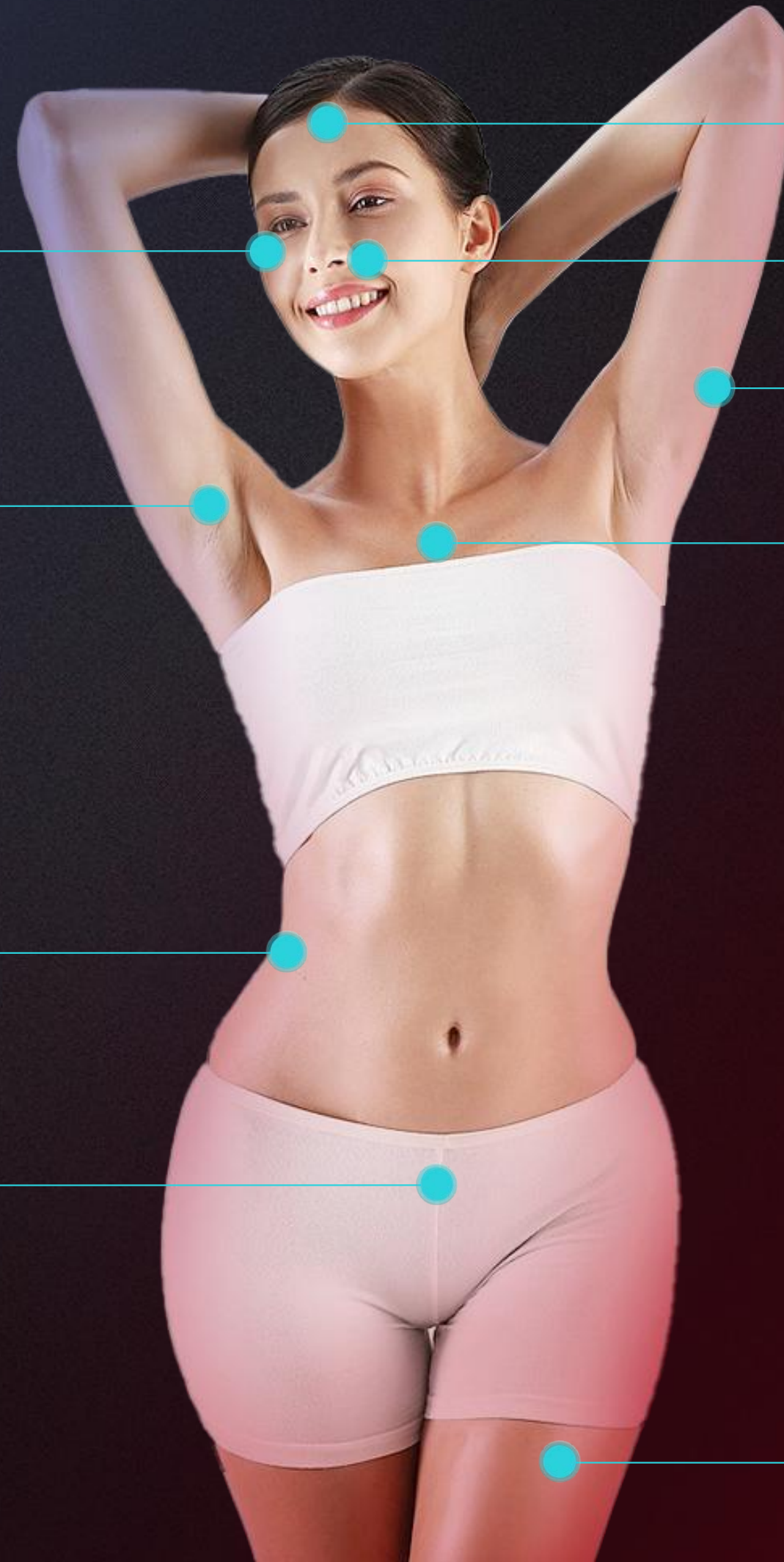

# *FAC Technology* *Fast Axis Collimator*

## **Collimator Lens Focus Beam**



Collimate the light path, reduce the divergence angle of the beam, and turn the divergent light path into a parallel light path. Smoother the edges of the laser beam.



# ***FAC Technology*** ***Fast Axis Collimator***



After addition of a collimator lens, beam focus

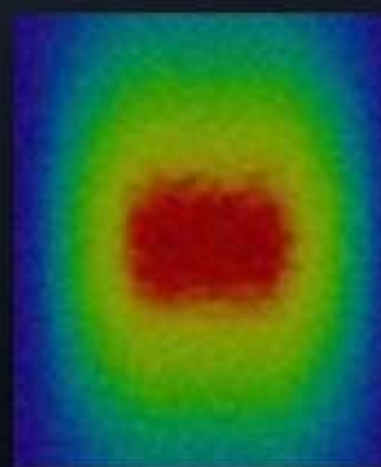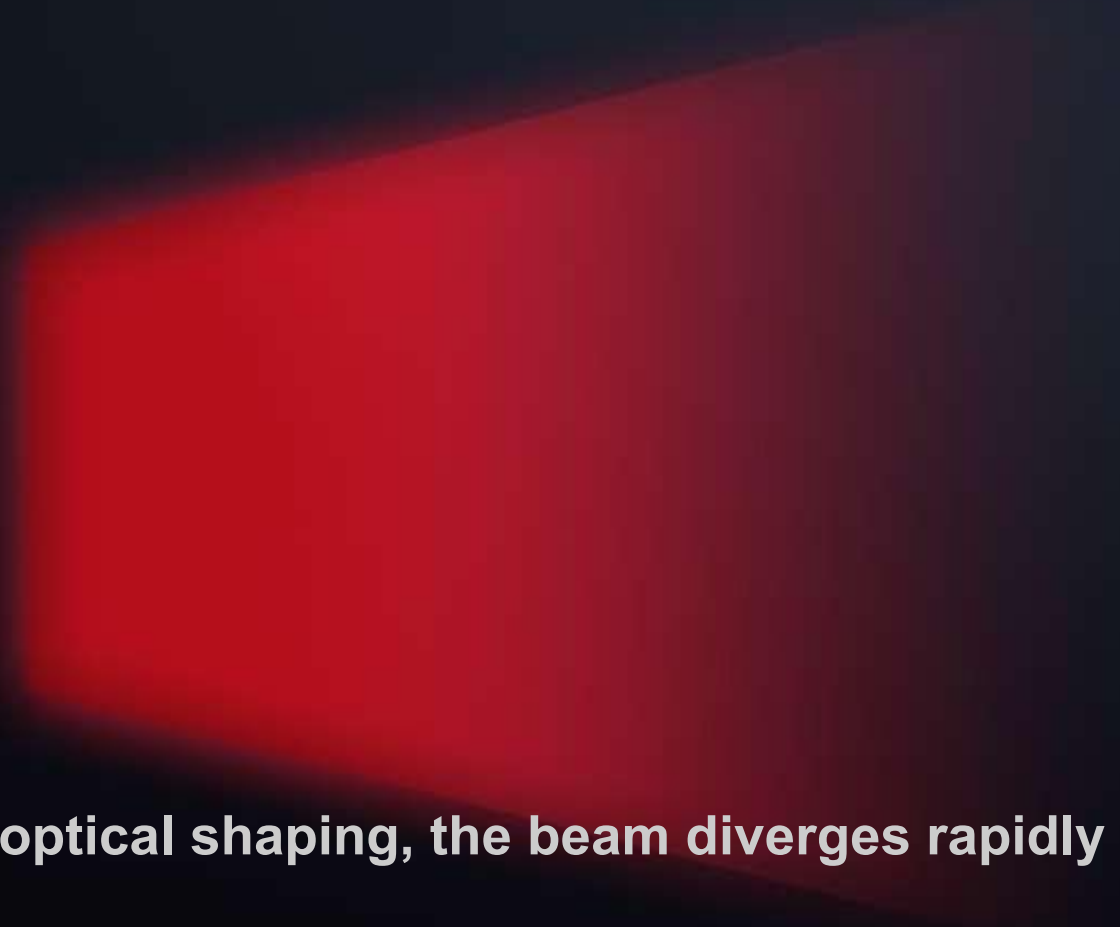

Without optical shaping, the beam diverges rapidly

***Large Spot***

**Super Large Spot Size:**

**5.2cm<sup>2</sup> 13.5\*38.5mm**

A unique and efficient new applicator with an especially large spot size of 5.2 cm<sup>2</sup>. This breakthrough applicator makes treatments much faster and pain free, significantly cutting treatment times and ensuring a more comfortable experience for patients.

**440% faster!**



# *Handle Screen*

## LCD Touch Screen



*Energy Adjustment*




*Pulse Width Adjustment*



*Frequency Adjustment*



*Real-time Temperature Monitoring*

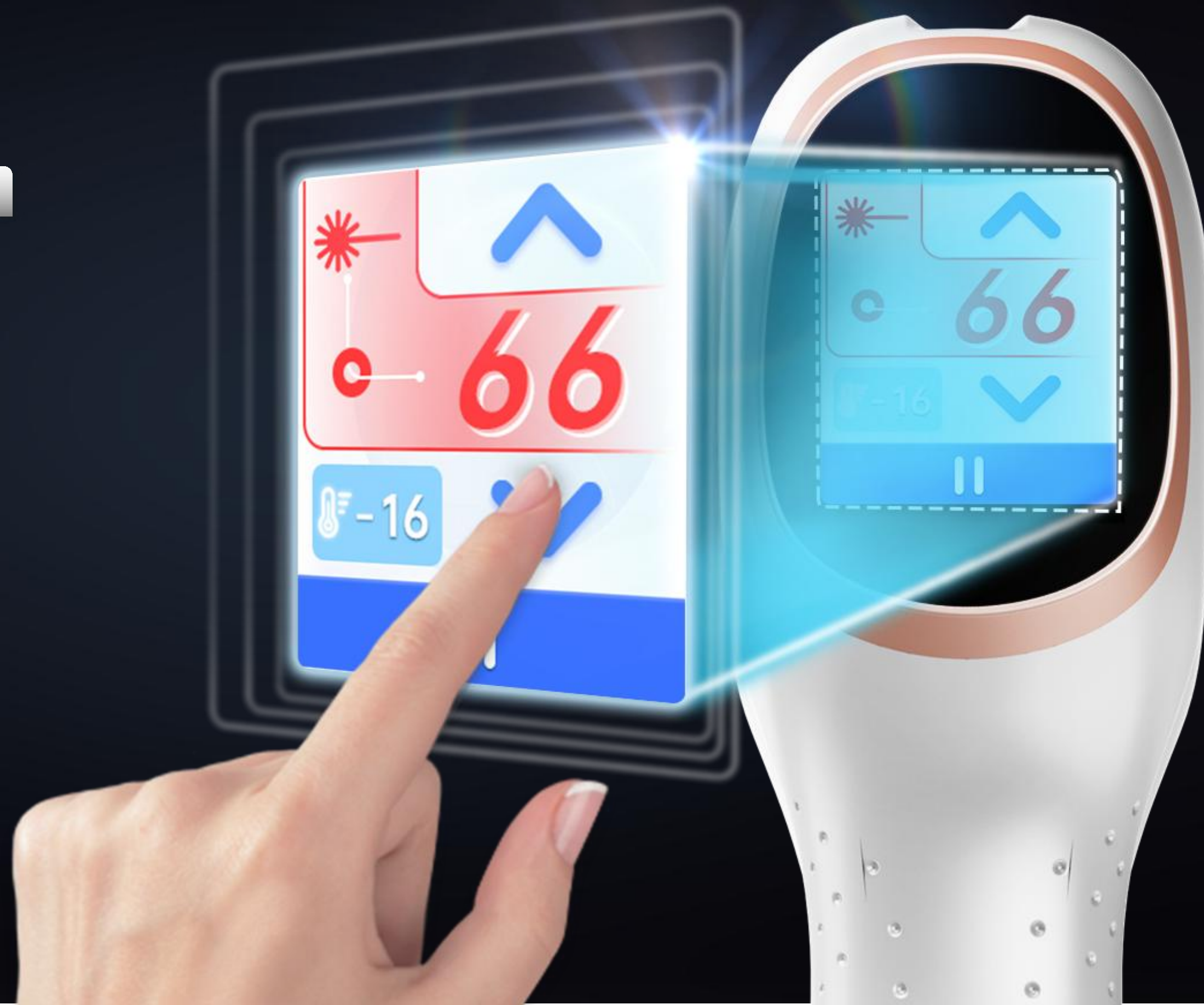

## ***Magnetic Heads***

# **Magnetic Heads, Strong And Safe.**

The magnetic connection is strong enough to resist most unintended disconnects, making sure the treatment is safe and perfect going.



***American Original***



## **American Original Import**

By Original American Coherent laser technology, the original bars are adopted in the laser generator to ensure that the output is stable and sufficient. The high quality bars guarantee a super zlong life span of the laser.

***Water Proof***

**Best water proof  
Results**



***Silent Wheel***

**High-End Hollow  
Silent Wheel**

*Stable  
Silent  
High Level*



**Low Noise**

## Say **"NO"** To Noise!

Silent working, make a more silent and convenient working environment for clients and operators.



## *Cooling Technology*

# **Handpieces With 360° Contact Cooling Strong Cooling, Painless Treatment.**

With strong TEC cooling technology, No pain during treatment, give patient a more comfortable treatment, and get a better treatment results.



# UI Interface



# Product Display

Handpiece Connector

Foot Switch

220V/110V Aviation Plug


Air Switch

Overflow

Water Input

Water Output

Water Filter Box

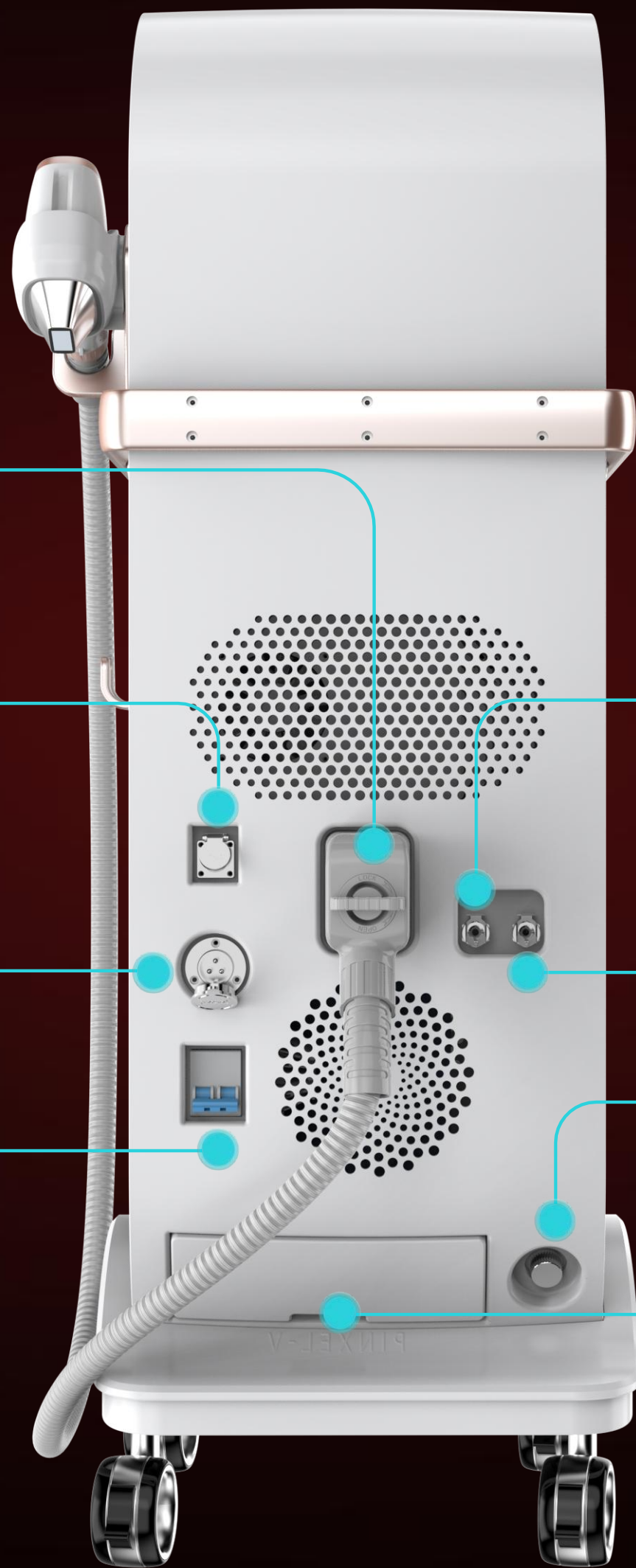

# Specifications

Wave Length:	<b>755nm 808nm 1064nm</b>
Fluence:	<b>5~80J/cm<sup>2</sup></b>
Pulse Width:	<b>5-100ms</b>
Spot Size:	<b>12.5*14.5mm 13.5*28.5mm 13.5*38.5mm Round 8.25mm</b>
Frequency:	<b>1~10Hz</b>
Handle Power:	<b>3000W</b>
Machine Power:	<b>6900W</b>
Cooling:	<b>TEC, air cooling + water cooling + semiconductor cooling</b>
Handpiece Cooling:	<b>-16°C</b>
Display:	<b>12 Inch multi color touch screen</b>



***Different Machine Color Available***



# Certificate

# US FDA EUROPE MEDICAL CE (TUV) APPROVED

The screenshot shows the FDA's 510(k) Premarket Notification database. The header includes the FDA logo and navigation links. The main content area displays the following information:

<b>Device Classification Name</b>	Powered Laser Surgical Instrument
<b>510(K) Number</b>	K191970
<b>Device Name</b>	Diode Laser Hair Removal System
<b>Applicant</b>	Beijing Sano Laser S&T Development Co., Ltd Room 7-201, No.1, Caids 3rd Street, Nancar Town, Shunyi District Beijing, CN
<b>Applicant Contact Correspondent</b>	Hongbo Zhang Microkn Business Consulting (Shanghai)Co., Ltd Room 1318, Block A, No 3698, Gonghexin Road, Jingan District Shanghai, CN 200435
<b>Correspondent Contact</b>	Huifang Zhao
<b>Regulation Number</b>	272.4010
<b>Classification Product Code</b>	3EX
<b>Date Received</b>	07/24/2019
<b>Decision Date</b>	04/06/2020
<b>Decision</b>	Substantially Equivalent (SESE)
<b>Regulation Medical Specialty</b>	General & Plastic Surgery
<b>510k Review Panel</b>	General & Plastic Surgery
<b>Type</b>	Traditional
<b>Reviewed By Third Party</b>	No
<b>Combination Product</b>	No



EBS medical Finland \* Korvenrannantie 25 A 30 \* 04320 Tuusula - Finland \* Y-tunnus: 3427432-5  
Phone: (+358) 45 3402360 / E-mail: infofinlanmd@ebs-medical.com / WEB: www.ebs-medical.fi

# Before & After



Before



After



Before



After



Before



After



Before



After

